



PROFINET

MBC3 Hook-up diagram ML Series instruments

PROFINET connection




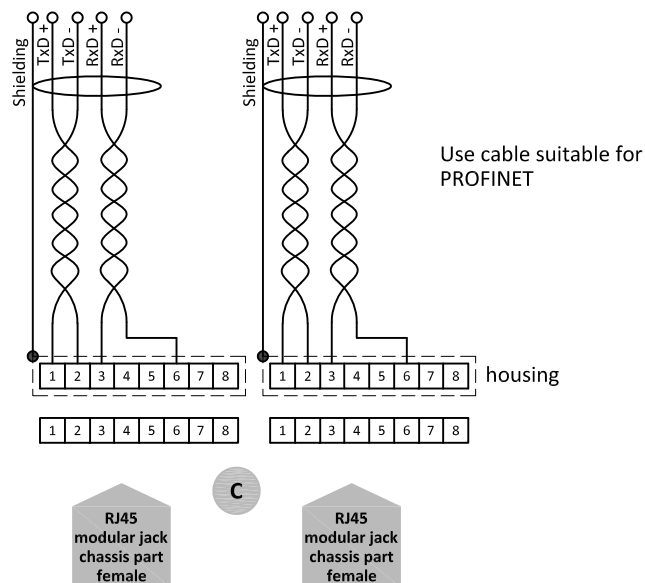


RS232 connection

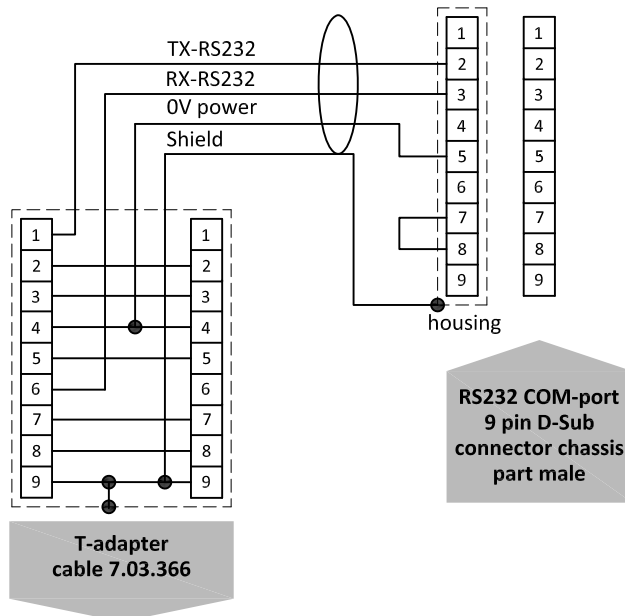

Model key explanation

For other explanation see 9.16.131

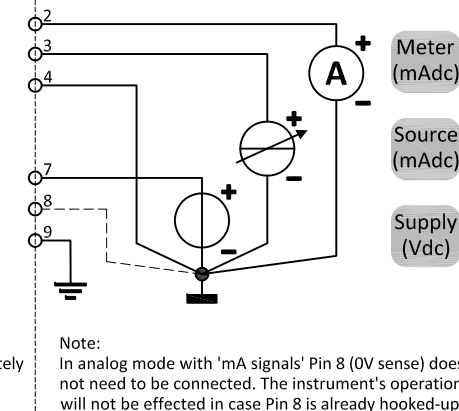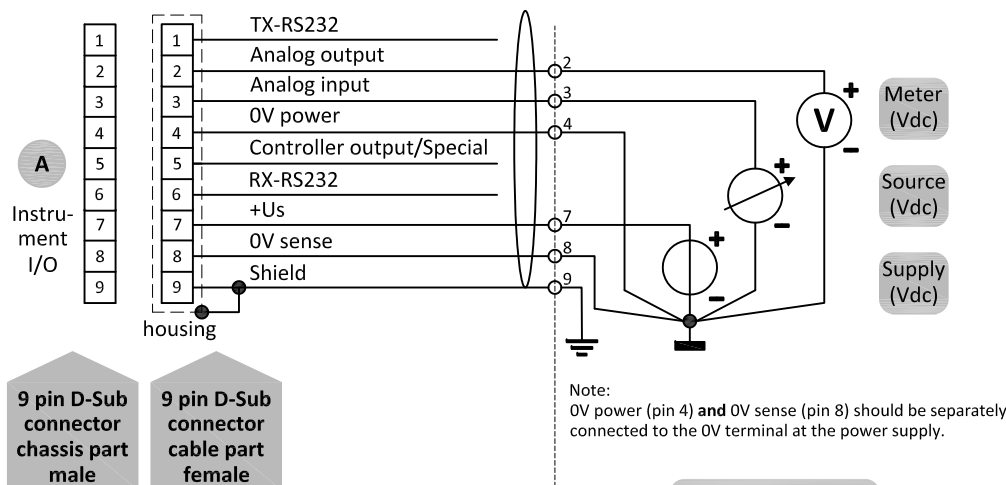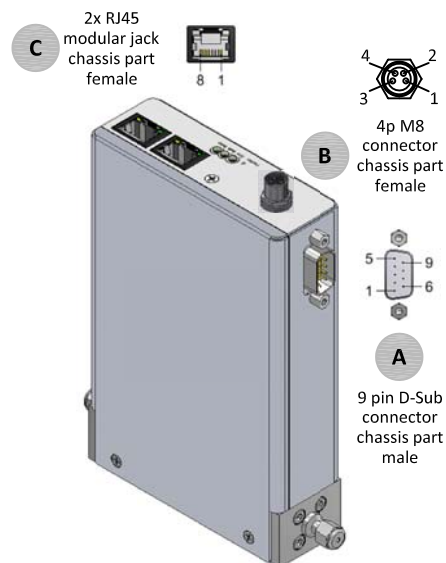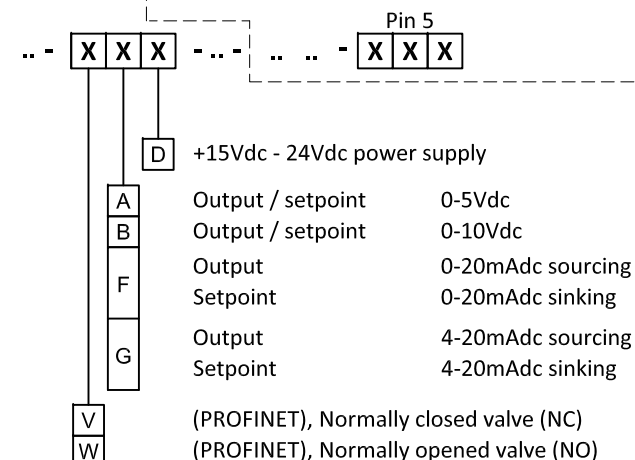

Note:
0V power (pin 4) and 0V sense (pin 8) should be separately connected to the 0V terminal at the power supply.

Note:
In analog mode with 'mA signals' Pin 8 (0V sense) does not need to be connected. The instrument's operation will not be effected in case Pin 8 is already hooked-up

**Analog operated
0-5 or 0-10Vdc**

**Analog operated
0-20 or 4-20mA dc**

Note:
For information about the optional B "Valve/Actuator" connector.
See doc.nr. 9.16.131 for more details.

Note:
Do not connect an external valve to instruments, set as meter.

Note:
When using a field bus or RS232, it is not possible to operate the instrument by using the setpoint signal of the analog D-sub connector without changing the value of parameter "control mode". See doc.nr. 9.17.023 for more details